
Hot-Cold Massage

Spot Therapy Vibrating Gel Pack

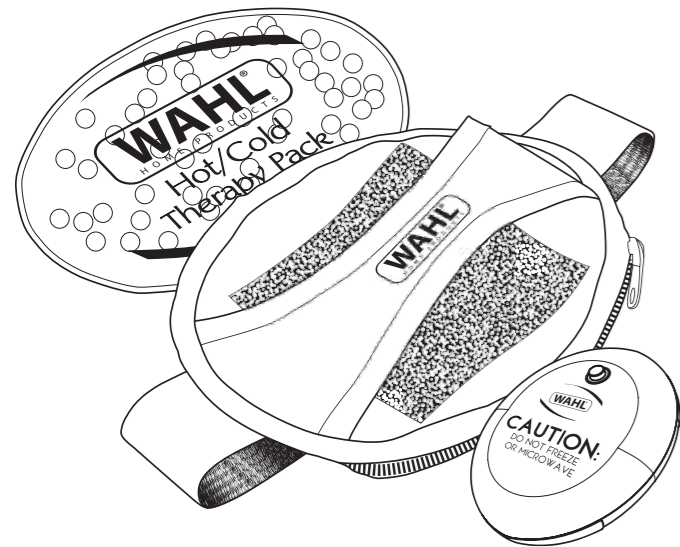The Brand Used By Professionals® Since 1919

© 2014 Wahl Clipper Corporation • Part No. 95810 • Printed in China
massagers.wahl.com

COLD THERAPY

Cold helps to reduce the swelling and pain caused by inflammation and relieves discomfort. May help with: temporary headache relief, fever reduction, toothache, sunburn, minor burns, sprains, minor cuts, beauty treatment, swelling reduction, scratches and abrasions. Contact your doctor if pain persists.

- Place gel pack in the freezer.
- After 1-2 hours, the gel pack will reach a therapeutic cold temperature. It is ideal to allow the gel pack to freeze overnight for optimal cold therapy.
- Cold will last for 20-30 minutes then remain cool.
- Do not apply the gel pack for longer than 20 minutes at one time.
- Wait at least 2 hours before beginning another cold therapy session.
- If the gel pack is uncomfortable during use, discontinue use immediately.

CERAMIC BEADS

Ceramic is a primary raw material commonly found in nature. These ceramic beads are non-toxic and are used in the gel pack to keep the pack hotter or colder longer. Do not ingest the ceramic beads. Follow the instructions in the cautions and warnings section if the gel pack becomes punctured or torn. If the ceramic beads are swallowed, contact a physician immediately.

OPERATING INSTRUCTIONS FOR MASSAGING DISC

1. Push button on the massaging disc to turn ON. Pushing the button a second time will turn the massaging disc OFF.
2. Place the massaging disc into the small zippered compartment of the hot-cold pack cover.
3. Apply the massaging hot-cold pack to the targeted area.

5

IMPORTANT SAFEGUARDS

When using an electrical appliance, basic precautions should always be followed, including the following:

READ ALL INSTRUCTIONS AND SAFEGUARDS BEFORE USING.

HOT-COLD GEL PACK

CAUTIONS AND WARNINGS

- Follow instructions carefully. Misuse can cause burns.
- Always be sure the gel pack is wrapped in cloth before applying to skin. We recommend using this gel pack with the included pouch to prevent injury. Direct application of this gel pack to the skin may cause burns or frostbite.
- Test the temperature by applying the gel pack to the inside of your wrist before using.
- Do not overheat. Overheating or excessive heating in the microwave may damage or destroy the gel pack. See chart for heating schedule.
- If any swelling of the gel pack is observed while heating, immediately turn off the microwave and wait for the swelling to go down before opening the microwave door.
- The gel pack and contents are safe when used as directed. If the gel pack is accidentally punctured or tears, discontinue using immediately. If the gel comes in contact with the skin, wash with soap and water. If skin becomes irritated or if the gel enters the eyes or mouth, contact a physician immediately.

TROUBLESHOOTING


If your massager fails to operate when switched "ON," check your batteries to see if they are properly installed. Positive end (center post) should be located as shown in battery housing. Make sure the switch is in the "ON" position. If it still does not operate, your batteries have probably exceeded their useful life. To maximize battery life, replace both AAA batteries at the same time.

BATTERY REPLACEMENT

When the massager's strength becomes weak, the battery needs to be replaced.

Do not mix old and new batteries. Do not mix alkaline, standard (carbon-zinc), or rechargeable (Ni-Cad, Ni-MH, etc.) batteries. Alkaline batteries made by any of the following vendors are preferred: Duracell™, Energizer™, or Rayovac™.

1. Press on the "key" and pry off the battery cover (Fig. 1).
2. Remove the batteries and insert two new AAA batteries. Make sure the right battery has the positive end up and the left battery has the positive end down (Fig. 2).
3. Replace the battery cover by inserting the bottom of the cover first then pressing on the upper part.

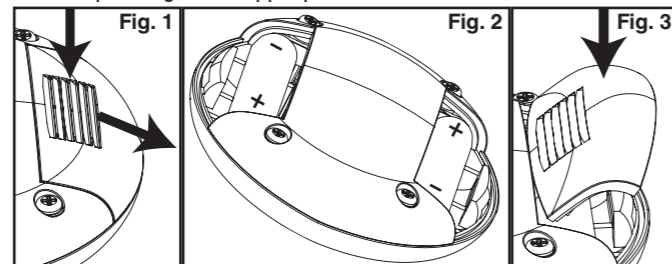6

- Check local laws and regulations for proper disposal of the gel pack.
- Be cautious when using the product on people whose skin may be more sensitive to temperature extremes, such as infants, young children or older persons.
- Persons with diabetes, circulatory problems, nerve damage, paralysis or sensitive skin should use only as directed by a physician.
- Consult your doctor before treating any serious injury or if pain persists.
- Do not use in conjunctions with other lotions, creams, or balms.
- For external use only.
- Use only as directed. Failure to follow these warnings may result in injury.
- Keep out of reach of children.

MASSAGING DISC

CAUTIONS AND WARNINGS

To reduce the risk of injury to persons:

1. **DO NOT MICROWAVE OR FREEZE THE MASSAGING DISC.** By doing so, you risk damaging both the massaging disc and your microwave. Placing batteries in a microwave oven may cause them to leak or explode and cause serious injury or property damage.
2. Do not use while bathing or in the shower.

7

CLEANING/STORAGE

IMPORTANT: If the massaging disc is going to be stored for a long period of time, be sure to remove the batteries.

To clean the massaging disc, simply wipe unit with a dry, soft cloth. Do not immerse in water or use any harsh cleaners or solvents.

To clean the protective cover, hand wash with a mild detergent. Lay flat and allow to dry completely before use.

OPERATING INSTRUCTIONS FOR HOT-COLD VIBRATING GEL PACK

1. Remove the gel pack from the microwave or freezer.
2. Open the velcro on the end of the protective cover. Slide the hot-cold gel pack into the protective cover. Close the velcro to secure the hot-cold gel pack in the protective cover.
3. Unzip the small zippered compartment in the middle of the gel pack cover.
4. Turn on the massaging disc.
5. Place the massaging disc into the small zippered compartment and close the zipper.
6. Apply the massaging hot-cold gel pack to the targeted area.

DO NOT MICROWAVE OR FREEZE THE MASSAGING DISC.

By doing so, you risk damaging both the massaging disc and your microwave.

Wahl® offers a full line of therapeutic massagers.


For more information or Product Registration, visit massagers.wahl.com

8

3. Do not place or store appliance where it can fall or be pulled into a tub or sink. Do not place in or drop into water or other liquid.
4. Close supervision is necessary when this appliance is used by, on, or near children or disabled persons.
5. Use this appliance only for its intended use, as described in this manual.
6. In case of pregnancy, diabetes or illness, consult your physician before using this massager.
7. Do not use massager on:
 - a. The front of the neck.
 - b. Genital areas.
 - c. Sleeping or unconscious persons.
8. Consult your doctor before using on swollen or inflamed areas.
9. Do not use this massager for more than 20 minutes at one time.

SAVE THESE INSTRUCTIONS For Household Use Only

USING HOT-COLD GEL PACK & MESSAGE TOGETHER MAY ADD EXTRA PAIN RELIEF

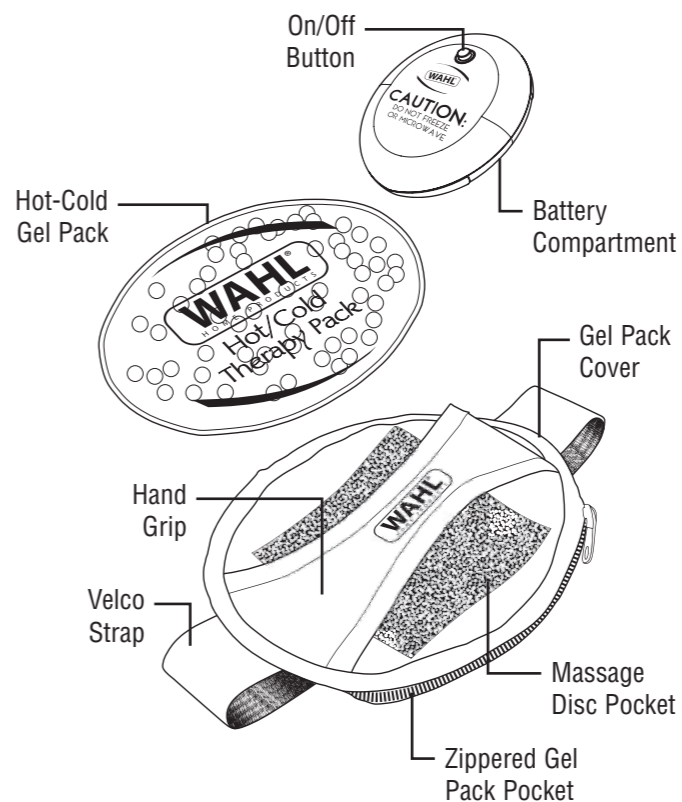8

OPERATING INSTRUCTIONS FOR HOT-COLD GEL PACK

HEAT THERAPY

Heat brings increased blood circulation to affected area. This helps to relax tight muscles and restore flexibility in a wide variety of ailments. May help with: temporary muscle ache and pain relief, neck and back pain, arthritis, muscle cramps, stomach ache, menstrual cramps, sinus pain, pelvic and leg pain. Contact your doctor if pain persists.

- Gel pack may be placed in microwave directly from freezer.

- Heating Instructions:

MICROWAVE OVEN HEATING SCHEDULE

Wattage	Heating Time
600	35 Seconds
800	30 Seconds
1000	20 Seconds
1300	15 Seconds

- Adjust heating times accordingly to suit the power of your microwave.

- Check the gel pack to make sure the desired temperature has been reached, and if needed, continue to heat in 10-second intervals.

IMPORTANT: Do not exceed heating times. To avoid injury do not exceed recommended heating times. In the event heating times are exceeded allow product to cool sufficiently before removing from microwave.

- Overheating will damage the product.

- Heat will generally last for 20-30 minutes then remain warm.

- Do not apply the product for longer than 20 minutes at one time.

- Wait at least 2 hours before beginning another heat therapy session.

- If the gel pack is uncomfortable during use, discontinue use immediately.

9

SERVICE AND REPAIR

This unit has no user-serviceable parts. If you have problems with your massager, **do not return it to store where purchased.** Send it prepaid with a note explaining the problem to:

Wahl Clipper Corporation
Attention: Repair Department
3001 Locust Street
Sterling, Illinois 61081
Telephone: (800) 767-9245

For products purchased outside the United States, contact your Wahl distributor.

2-YEAR LIMITED WARRANTY

Valid in U.S.A. Only

If your Wahl product fails to operate satisfactorily within two (2) years from the date of original purchase or receipt as a gift, Wahl will, at our option, repair or replace it and return it to you FREE OF CHARGE. Please mail your product prepaid to Wahl Clipper Corporation, 3001 Locust Street, Sterling, Illinois 61081, Attention: Repair Department, with a note describing the defect and the date of original purchase or receipt as a gift (if in warranty). Products out of warranty will be repaired at our standard repair charges.

This warranty is void and no repair or replacement will be made under this warranty or otherwise if:

1. The product is modified in any manner or repaired by anyone other than Wahl Clipper Corporation or Wahl authorized service center.
2. The product has been used commercially or subjected to unreasonable use.

This warranty gives you specific legal rights and you may also have other rights which may vary from state to state. No other written express warranty is given. Any implied warranty, including any warranty of merchantability which may arise from purchase or use is limited to the years period provided in the express warranty. Some states may not allow such limitation, so it may not apply to you. No responsibility is assumed for incidental or consequential damages of any type or for the use of unauthorized attachments. Some states may not allow the exclusion or limitation of incidental or consequential damages, so the limitation or exclusion may not apply to you.

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WAHL GmbH