
Neck & Back Massage

Deep-Vibration Heated Massage Wrap

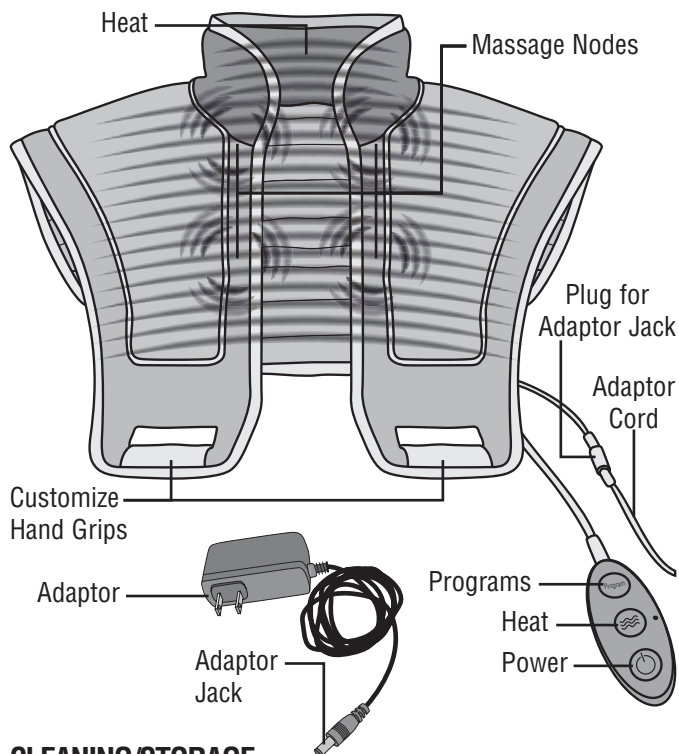CLEANING/STORAGE

Always unplug massager before cleaning.

To clean, simply wipe unit with a dry, soft cloth. Do not immerse in water or use any harsh cleaners or solvents. When storing unit, make sure cord is not twisted or kinked.

For more information or Product Registration, visit massagers.wahl.com

1

OPERATING INSTRUCTIONS

Wahl's Neck & Back Massage Wrap is specifically designed to fit the neck, shoulder, and upper back areas where users have the most physical signs of pain and tension. Drape the Neck & Back Massage Wrap comfortably on your shoulders making sure the neck flap is properly in place. This unit is weighted to help stay in place. The hand grips help to customize the massage experience.

1. Insert adaptor jack into the plug on the left side of the unit.
2. Plug into a AC 120-volt outlet.
3. To turn the unit on, press the power button on the controller once. The unit will not turn on until the power button is pushed.
4. To turn the heat on, press the heat button once. Press again to turn it off. The Wahl Neck & Back Massage Wrap has been designed to provide optimum heat.
5. To change the programs, press the program button. Select from multiple settings by pressing the programs button successively until desired program is found.
6. To turn the unit off, hold the power button for 3 seconds. This unit automatically turns off after 20 minutes.

SERVICE AND REPAIR

This unit has no user-serviceable parts. If you have problems with your massager, **do not return it to store where purchased.** Send it prepaid with a note explaining the problem to:

Wahl Clipper Corporation
Attention: Repair Department
3001 Locust Street
Sterling, Illinois 61081
Telephone: (800) 767-9245

For products purchased outside the United States, contact your Wahl distributor.

2

IMPORTANT SAFEGUARDS

When using an electrical appliance, basic precautions should always be followed, including the following:

READ ALL INSTRUCTIONS AND SAFEGUARDS BEFORE USING.

DANGER

To reduce the risk of electric shock:

1. Do not reach for an appliance that has fallen into water. Unplug immediately.
2. Do not use while bathing or in the shower.
3. Do not place or store appliance where it can fall or be pulled into a tub or sink. Do not place in or drop into water or other liquid.
4. Always unplug this appliance from the electrical outlet immediately after using and before cleaning.
5. Never use pins or other metallic fasteners with this appliance.
6. Carefully examine the covering before each use. Discard the appliance if the covering shows any sign of deterioration, such as checking, blistering, or cracking.
7. Keep Dry – Do not operate in a wet or moist condition.

WARNING

To reduce the risk of burns, fire, electric shock or injury to persons:

1. This appliance should never be left unattended when plugged in. Unplug from outlet when not in use, and before putting on or taking off parts.
2. Do not operate this appliance under a blanket or pillow. Excessive heating can occur and cause fire, electric shock or injury to persons.
3. Close supervision is necessary when this appliance is used

3

- by, on, or near children or disabled persons.
4. This massager contains a motor. To avoid risk of injury, it must not be used within six inches of an implantable cardiac device.
5. Use heated surfaces carefully. May cause serious burns. Do not use over insensitive skin areas or in the presence of poor circulation. The unattended use of heat by children or incapacitated persons may be dangerous.
6. Use this appliance only for its intended use, as described in this manual.
7. Never operate this appliance if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Return to Wahl for examination and repair.
8. Do not carry this appliance by supply cord or use cord as a handle.
9. Keep the cord away from heated surfaces.
10. Never drop or insert any object into any opening on this appliance.
11. Do not use outdoors.
12. Do not operate where aerosol (spray) products are being used or where oxygen is being administered.
13. In case of pregnancy, diabetes or illness, consult your physician before using this massage wrap.
14. Do not use massage wrap on:
 - a. Sensitive skin, swollen or inflamed areas, in the presence of poor circulation, where skin eruptions are present.
 - b. The front of the neck.
 - c. Genital areas.
 - d. Sleeping or unconscious persons.
15. Do not use this massage wrap for more than 20 minutes at one time.

4

16. Before you unplug from outlet, hold the power button for 3 seconds to turn OFF and press the heat button to turn OFF, then remove plug from outlet.
17. Do not use on an infant or invalid or on a sleeping or unconscious person. Do not use on insensitive skin or on a person with poor blood circulation. Check the skin in contact with the heated area of the appliance frequently to reduce the risk of blistering.
18. When storing this product, do not crush and avoid sharp folds.

SAVE THESE INSTRUCTIONS For Household Use Only

1-YEAR LIMITED WARRANTY

Valid in U.S.A. Only

If your Wahl product fails to operate satisfactorily within one (1) year from the date of original purchase or receipt as a gift, Wahl will, at our option, repair or replace it and return it to you FREE OF CHARGE. Please mail your product prepaid to Wahl Clipper Corporation, 3001 Locust Street, Sterling, Illinois 61081, Attention: Repair Department, with a note describing the defect and the date of original purchase or receipt as a gift (if in warranty). Products out of warranty will be repaired at our standard repair charges.

This warranty is void and no repair or replacement will be made under this warranty or otherwise if:

1. The product is modified in any manner or repaired by anyone other than Wahl Clipper Corporation or Wahl authorized service center.
2. The product has been used commercially or subjected to unreasonable use.

This warranty gives you specific legal rights and you may also have other rights which may vary from state to state. No other written express warranty is given. Any implied warranty, including any warranty of merchantability which may arise from purchase or use is limited to the years period provided in the express warranty. Some states may not allow such limitation, so it may not apply to you. No responsibility is assumed for incidental or consequential damages of any type or for the use of unauthorized attachments. Some states may not allow the exclusion or limitation of incidental or consequential damages, so the limitation or exclusion may not apply to you.

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Wahl Clipper Corporation, Attention: Repair Dept.